


Reversed Phase SPE

VertiPak™ C4

- VertiPak™ C4 SPE is moderately non-polar, hydrophobic sorbent, manufactured using a trifunctional silane and endcapped to reduce the polar secondary interactions associated with surface silanol groups
- Less retentive than VertiPak™ C8 based on non-polar interaction and slightly more polar than C18 or C8 due to short chain length of C4
- Typical applications include compounds which are too strongly retained on C18 or C8 e.g. analgetics in blood
- Tube format is available in 7 tube sizes
- Cartridge format is available in 300mg, 600mg and 900mg size
- Tubes or cartridges are packaged in zip sealing bag protective from moisture and light
- Includes Certificate of Analysis



Chemical structure

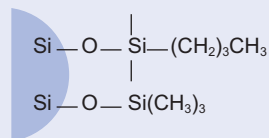

Chemical structure of C4 (Butyl) silane and trimethyl silyl group, covalently bonded to the surface of a silica particle.

Specifications

Phase	C4 (Butyl)
Base	Irregular-shaped Silica
Particle Size (µm)	50
Pore Size (Å)	60
Surface Area (m ² /g)	500
Carbon Load (%)	8
Bonding	Trifunctional
Endcap	Yes
pH Stability	2-9

Ordering Information

Bed Weight	Tube Size	QTY	Part No.
VertiPak™ C4 SPE Tubes			
50mg	1mL	100	0128-0236
50mg	3mL	50	0128-0245
100mg	1mL	100	0128-0336
100mg	3mL	50	0128-0345
100mg	10mL	30	0128-0364
200mg	3mL	50	0128-0445
200mg	10mL	30	0128-0464
500mg	3mL	50	0128-0645
500mg	6mL	30	0128-0654
500mg	10mL	30	0128-0664
1,000mg	6mL	30	0128-0954
2,000mg	12mL	30	0128-1074
5,000mg	25mL	20	0128-12A2
5,000mg	35mL	20	0128-1282
10,000mg	60mL	16	0128-1391
VertiPak™ C4 SPE Cartridges			
300mg		25	0128-0503
300mg		50	0128-0505
300mg		100	0128-0506
600mg		25	0128-0713
600mg		50	0128-0715
600mg		100	0128-0716
900mg		25	0128-0823
900mg		50	0128-0825
900mg		100	0128-0826